


SRM 870 Test Mix on Viva C18

1. uracil
2. toluene
3. ethylbenzene
4. quinizarin
5. amitriptyline

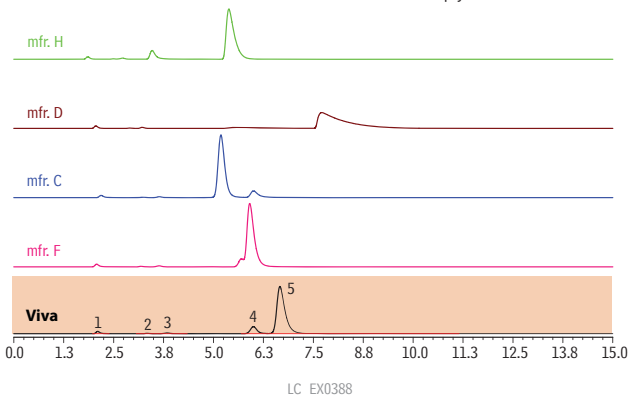

Sample:

Inj.: 10 μ L HPLC Performance (NIST 870) Test Mix (cat.# 31699)

Column: Viva C18

Cat.#: 9514562
Dimensions: 150 x 2.1mm wide pore
Particle size: 5 μ m
Pore size: 300 \AA

Conditions:

Mobile phase: water:methanol
(25:75 v/v)
Flow: 0.2mL/min.
Det.: UV @ 254nm